

Hi-Boy Whirlpools

Use



Division of Acorn Engineering Company®

H Series 60 Gallon Stationary Whirlpool

H-60-S



Whitehall H-60-S Stationary Whirlpool is recommended for shallow or deep body immersion offering complete coverage of the lower extremities as well as the lower back and hip.

Tank Capacity is 60 gallons.

Whirlpool Tank is 36"L x 20"W x 28"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup. Tank is also provided with an extended stainless steel skirt flush to the floor line to conceal internal drain piping. A 2" Combination Drain and Overflow, and a Filler Spout are provided.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 60-gallon capacity, 36"L x 20"W x 28"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a stainless steel skirt flush to the floor line. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Provide a 2" combination drain and overflow, and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit **www.whitehallmfg.com** for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/03/18




Hi-Boy Whirlpools

MODEL NUMBER AND OPTION SELECTIONS Series H Hi-Boy Series **Capacity Type** ✓ -S Stationary Whirlpool **Turbine Electrical Requirements** ☐ 115VAC / 50Hz ☐ 115VAC / 60HZ ☐ 230VAC / 50Hz ☐ 230VAC / 60HZ Country of Use

Refer to supplemental MXT15 sheet for details.

For Mixing Valve:

Suffix O	ptions (See Whirlpool Options Section for Details)
-AHC1	Adjustable High Chair with Casters
-AHC2	Adjustable High Chair with Rubber Tips
-ARS	Arm Rest Support
🔲 -AS1	Adjustable Suspension Seat
-CTS1	Combination Table with Side Seats
🔲 -EG	Enviro Glaze Color Powder Coating -
	Specify Color:
	(Special Finishes Consult Factory)
🔲 -GFI	Ground Fault Interrupter
-PCS	Pump Strainer Cover
🔲 -RHG	Rubber Hand Grips
🔲 -SDO	Separate Drain and Overflow Assembly
☐ -TTS1	Tank Top Seat

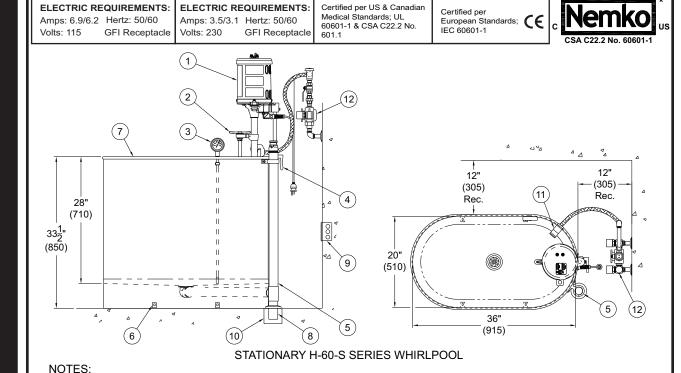

3. Thermometer 4. Turbine Raising & Lowering Device

1. 1/2" HP Jet Pump Motor

2. Variable Pressure Control

Overflow Assembly

- "L" Brackets (4) Four Each
- 7. Stainless Steel Tank
- 8. 2" O.D. x 5" Tailpiece

5. Combination Drain &

- 9. 3 Pole GFI Receptacle by Others
- 10. Floor Drain by Others

MARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

- 11. Filler Spout
- 12. Mixing Valve, Shown For Reference Only

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

			
SELECTION SUMMARY & APPROVAL FOR MANUFACTURING			
Model Number & Options		Quantity	
Company		Date	
Contact			
Approval for Manufacturing/Signature			

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862